



Simec AG

Areal Bleiche West Postfach 413 CH - 4800 Zofingen
Telefon +41 62 752 83 08 Fax +41 62 752 83 09
E-Mail: info@simec.ch www.simec.ch

Test report

Number **C 34811**

Object of analysis: Oxymetholone

Customer: Lightbest International Limited

Samples: 1

Batch: -

Receipt: 05.08.2015

Subject: Assay, Identification

Method(s): HPLC

Sample / Analysis:	Method	Status	Result	Specification	complies
<u>Oxymetholone:</u>					
Identification	HPLC	N	complies	complies	yes
Content	HPLC	N	101.6 %	$\geq 98.5 \%$	-

Method-status: G = GMP A = accredited V = generally validated P = validated on product N = not validated
E = external lab

Remarks:

Signatures: controlled & released:

Staff:

completed:


06.10.2015

page 1 of 1

Report steroid analysis

HPLC: Nr. 6
 Analysendatum: 05.10.2015 16:39:41
 Vial: 1
 Probenidentifikation: Blank
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Säule Nr. 85, 100-5 C18 250 x 4.6 mm
 Detektion: UV 240nm

Druckdatum: 06.10.2015 11:41:57
 Injektvolumen: 20

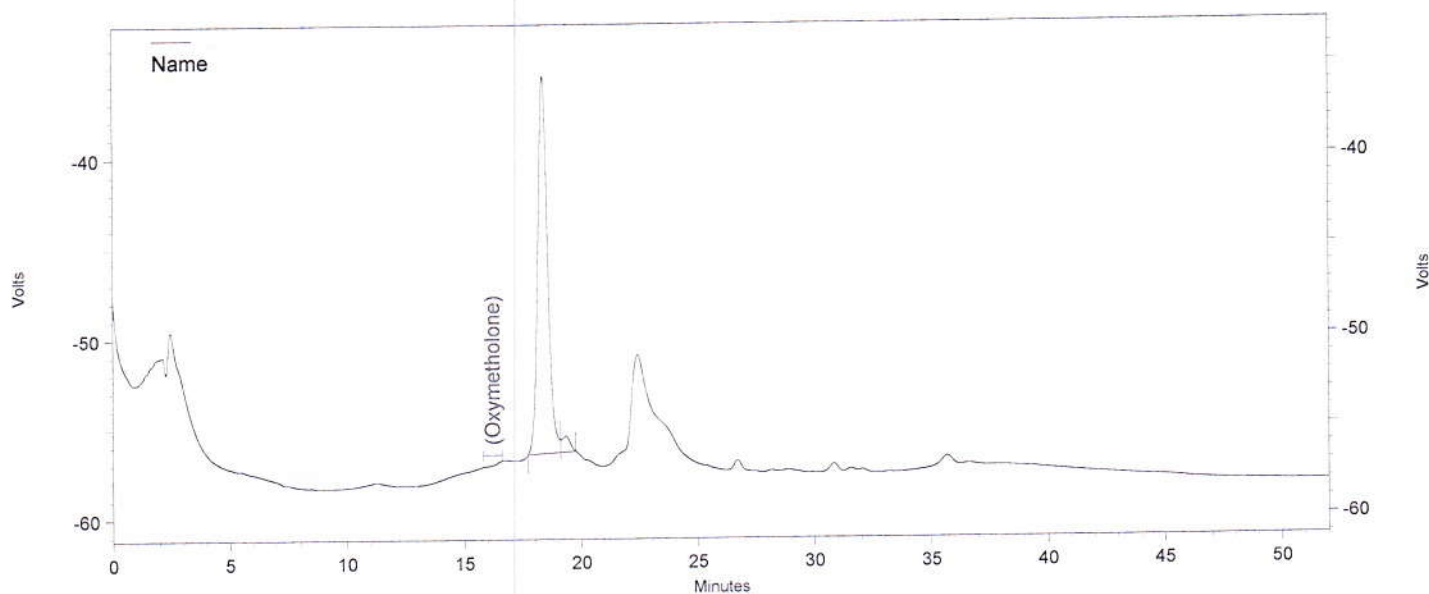

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone			0.000 BDL	
	18.387	2423745	0.000	96.645
	19.347	84150	0.000	3.355
Totals		2507895	0.000	100.000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
 Steroide\05.10.15\05.10.2015 16-38-06_Blank.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 05.10.2015 20:13:16
Vial: 2
Probenidentifikation: STD Oxymetholone
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 13:21:05
Injektvolumen: 20

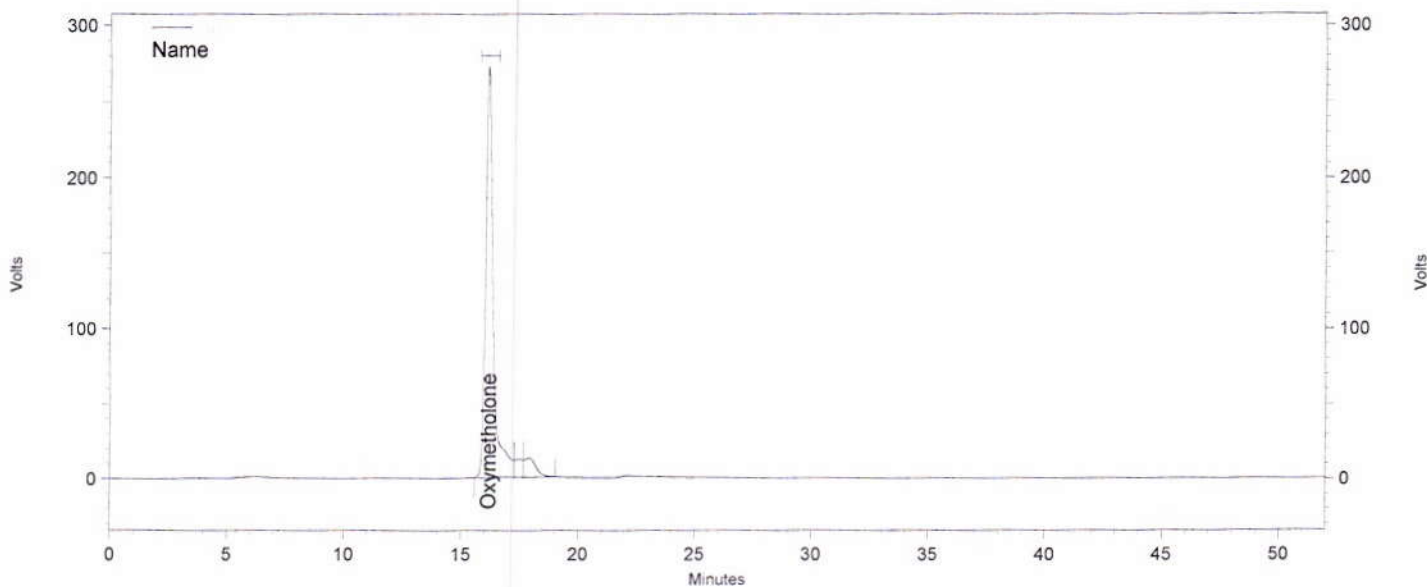

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	16.167	24727431	0.000	89.548
	17.483	1139987	0.000	4.128
	17.930	1746115	0.000	6.323
Totals		27613533	0.000	100.000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\05.10.2015 20-12-02_STD Trenbolone Acetat-Rep4.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 05.10.2015 17:33:02
Vial: 2
Probenidentifikation: STD Oxymetholone
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:43:28
Injektvolumen: 20

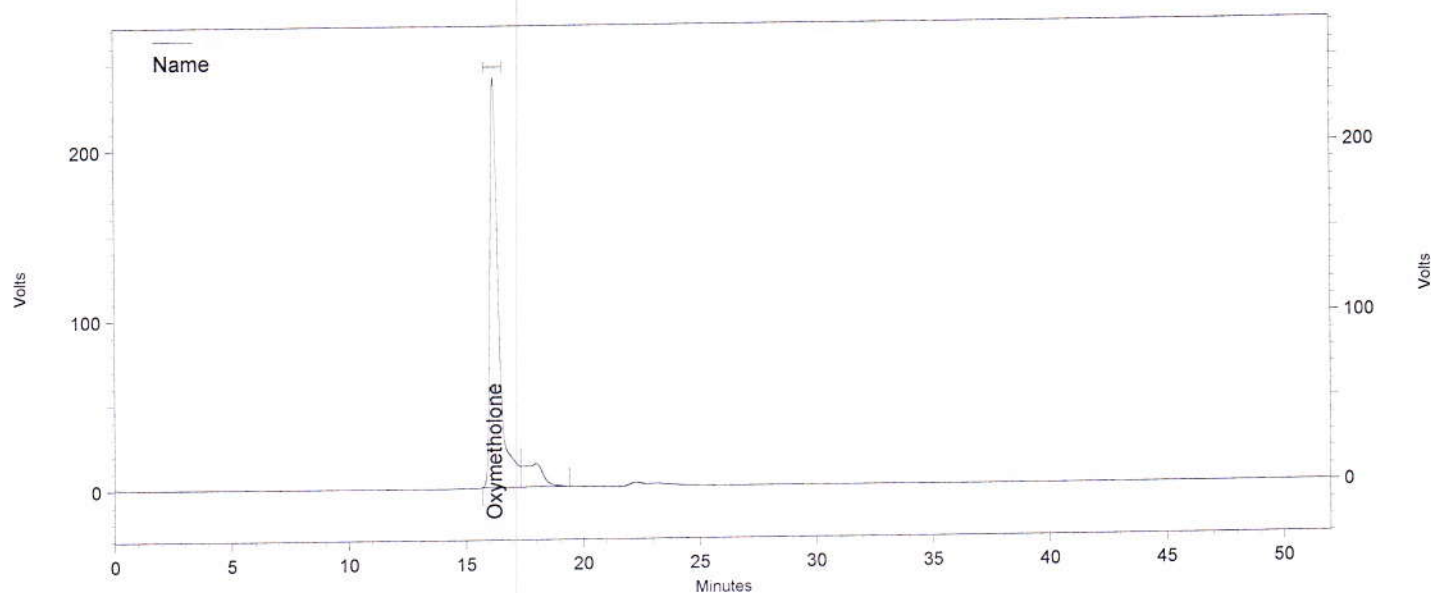

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	16.203	24562780	0.000	88.498
	18.000	3192479	0.000	11.502


Totals		27755259	0.000	100.000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\05.10.2015 17-31-48 STD Trenbolone Acetat-Rep1.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 05.10.2015 18:26:26
Vial: 2
Probenidentifikation: STD Oxymetholone
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:43:45
Injektvolumen: 20

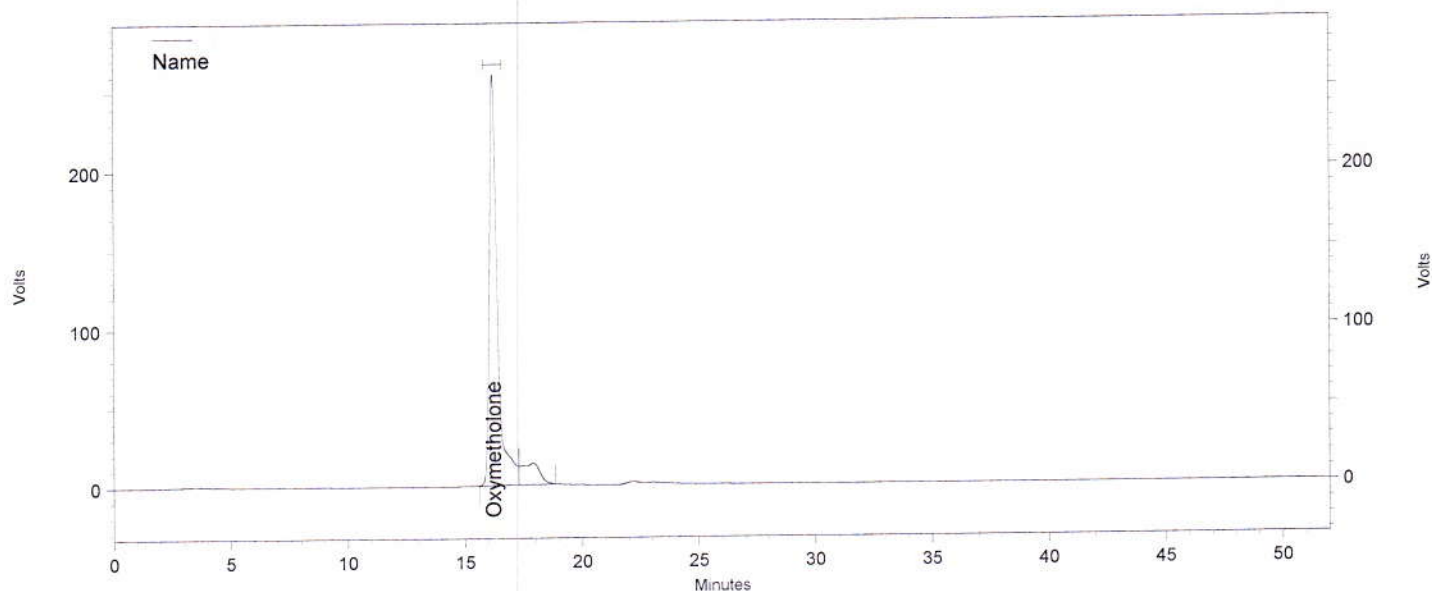

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	16.200	24729868	0.000	89.322
	17.920	2956297	0.000	10.678


Totals		27686165	0.000	100.000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\05.10.2015 18-25-13_STD Trenbolone Acetat-Rep2.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
 Analysendatum: 05.10.2015 19:19:49
 Vial: 2
 Probenidentifikation: STD Oxymetholone
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Säule Nr. 88, 100-5 C18 250 x 4.6 mm
 Detektion: UV 240nm

Druckdatum: 06.10.2015 11:44:05
 Injektvolumen: 20

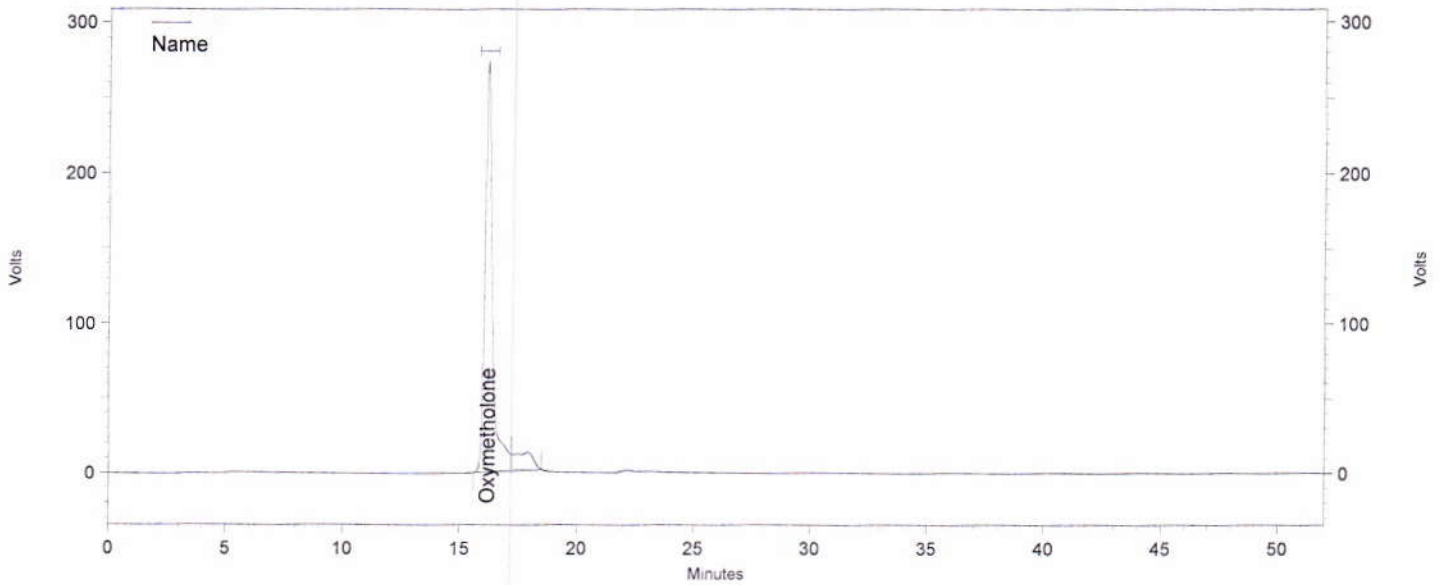

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	16.177	24394523	0.000	90.082
	17.927	2685969	0.000	9.918
Totals		27080492	0.000	100.000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
 Steroide\05.10.15\05.10.2015 19-18-33_STD Trenbolone Acetat-Rep3.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
 Analysendatum: 05.10.2015 21:59:58
 Vial: 2
 Probenidentifikation: STD Oxymetholone
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Säule Nr. 88, 100-5 C18 250 x 4.6 mm
 Detektion: UV 240nm

Druckdatum: 06.10.2015 11:44:14
 Injektvolumen: 20

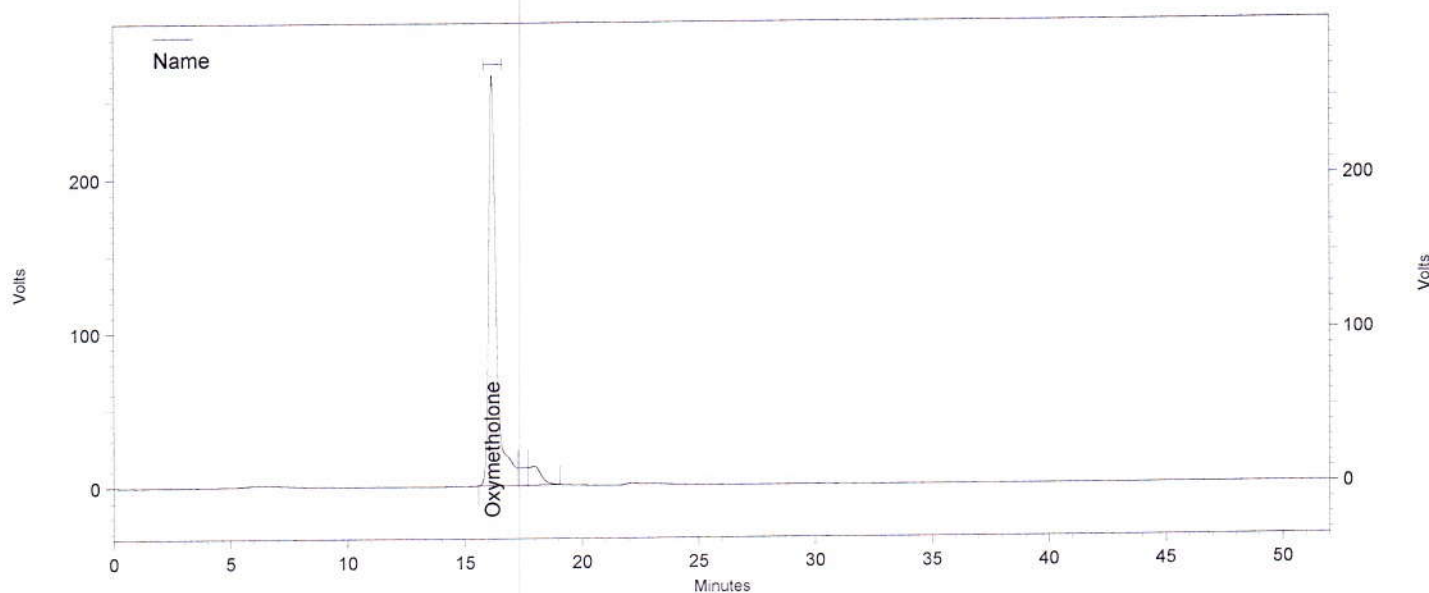

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	16.157	24761554	0.000	89.624
	17.517	1138368	0.000	4.120
	17.953	1728296	0.000	6.256
Totals		27628218	0.000	100.000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
 Steroide\05.10.15\05.10.2015 21-58-43 STD Trenbolone Acetat-Rep6.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
 Analysendatum: 05.10.2015 21:06:36
 Vial: 2
 Probenidentifikation: STD Oxymetholone
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Säule Nr. 88, 100-5 C18 250 x 4.6 mm
 Detektion: UV 240nm

Druckdatum: 06.10.2015 11:44:19
 Injektvolumen: 20

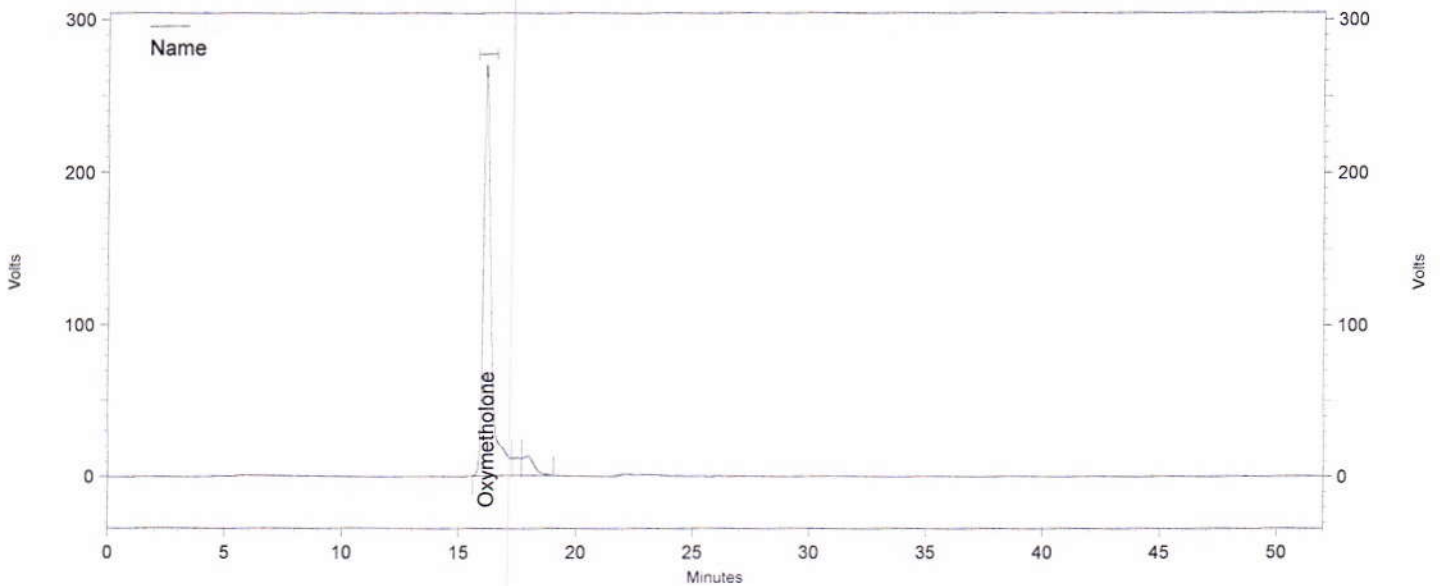

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	16.163	24821472	0.000	89.554
	17.497	1137338	0.000	4.103
	17.950	1758073	0.000	6.343
Totals		27716883	0.000	100.000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
 Steroide\05.10.15\05.10.2015 21-05-22_STD Trenbolone Acetat-Rep5.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 05.10.2015 22:53:19
Vial: 1
Probenidentifikation: Blank
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:44:26
Injektvolumen: 20

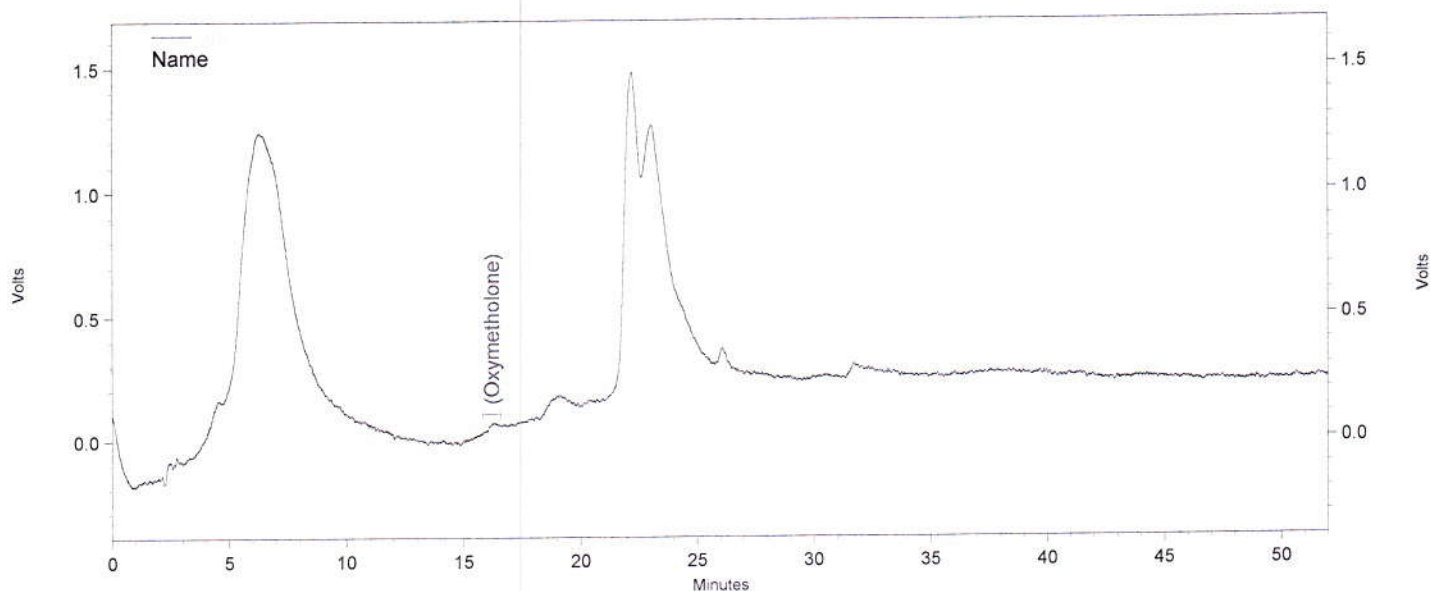

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone			0.000 BDL	
Totals				

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\05.10.2015 22-52-03_Blank.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 05.10.2015 23:46:44
Vial: 3
Probenidentifikation: C34811 A
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:44:44
Injektvolumen: 20

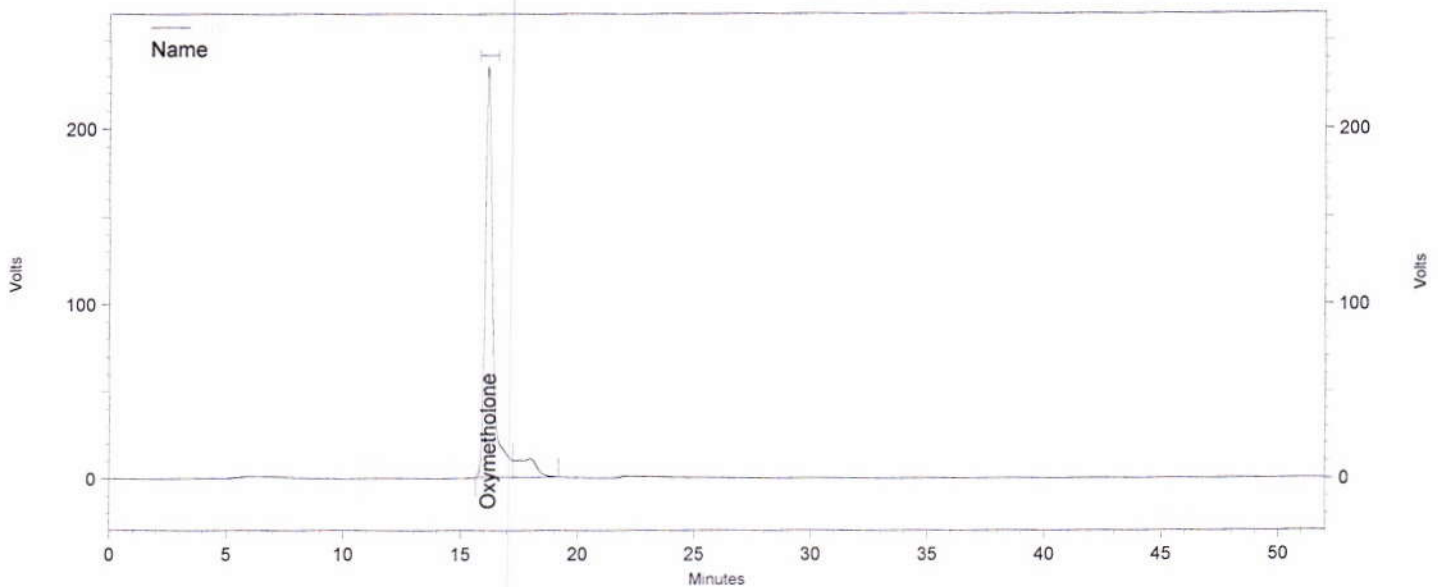

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	16.173	21612134	0.000	89.289
	17.993	2592584	0.000	10.711


Totals		24204718	0.000	100.000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\05.10.2015 23-45-32_C35895 A-Rep1.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 06.10.2015 00:40:10
Vial: 3
Probenidentifikation: C34811 A
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:44:59
Injektvolumen: 20

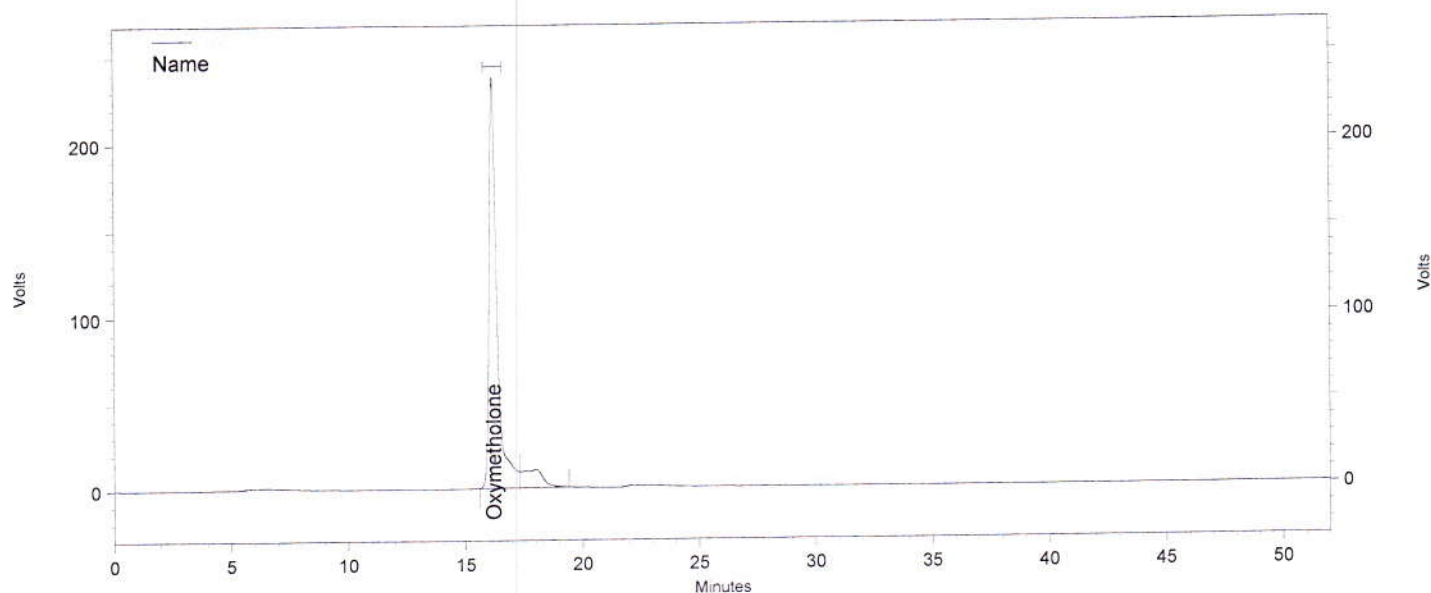

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	16.180	21935185	0.000	89.601
	18.007	2545757	0.000	10.399


Totals		24480942	0.000	100.000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\06.10.2015 00-38-57_C35895 A-Rep2.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 06.10.2015 01:33:35
Vial: 3
Probenidentifikation: C34811 A
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:45:18
Injektvolumen: 20

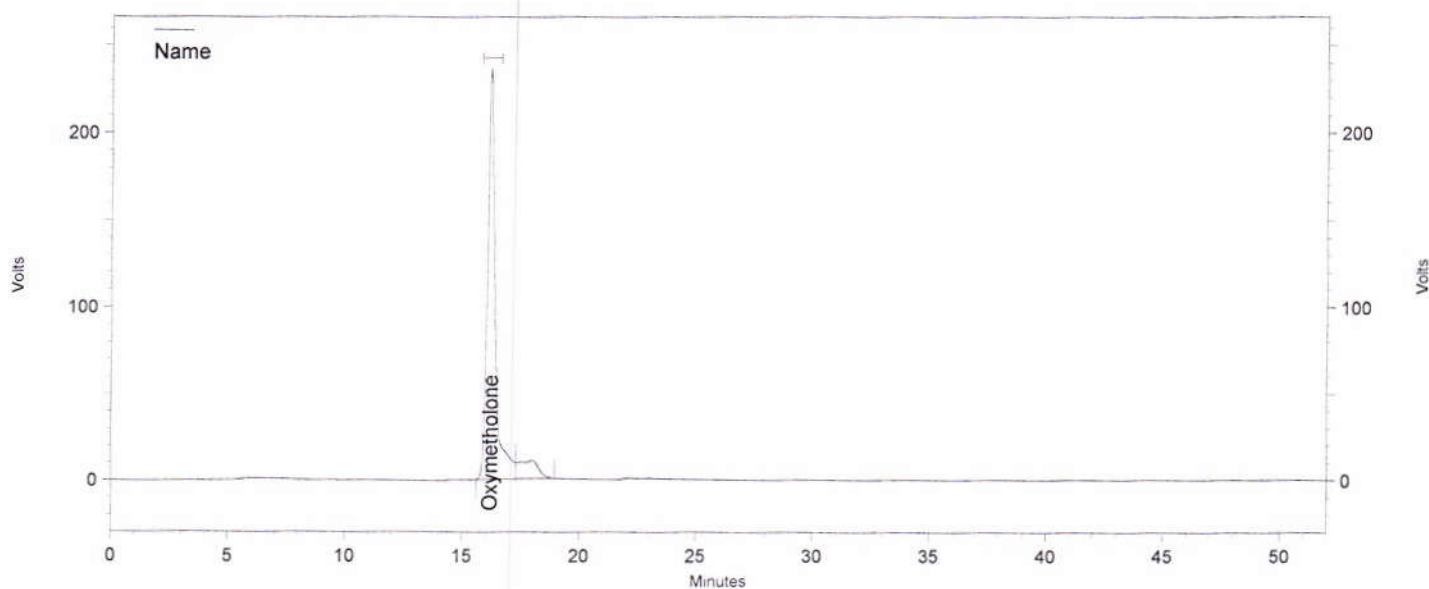

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	16.180	21771145	0.000	89.691
	18.007	2502227	0.000	10.309
Totals		24273372	0.000	100.000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\06.10.2015 01-32-23_C35895 A-Rep3.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 06.10.2015 02:26:57
Vial: 4
Probenidentifikation: C34811 B
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:46:08
Injektvolumen: 20

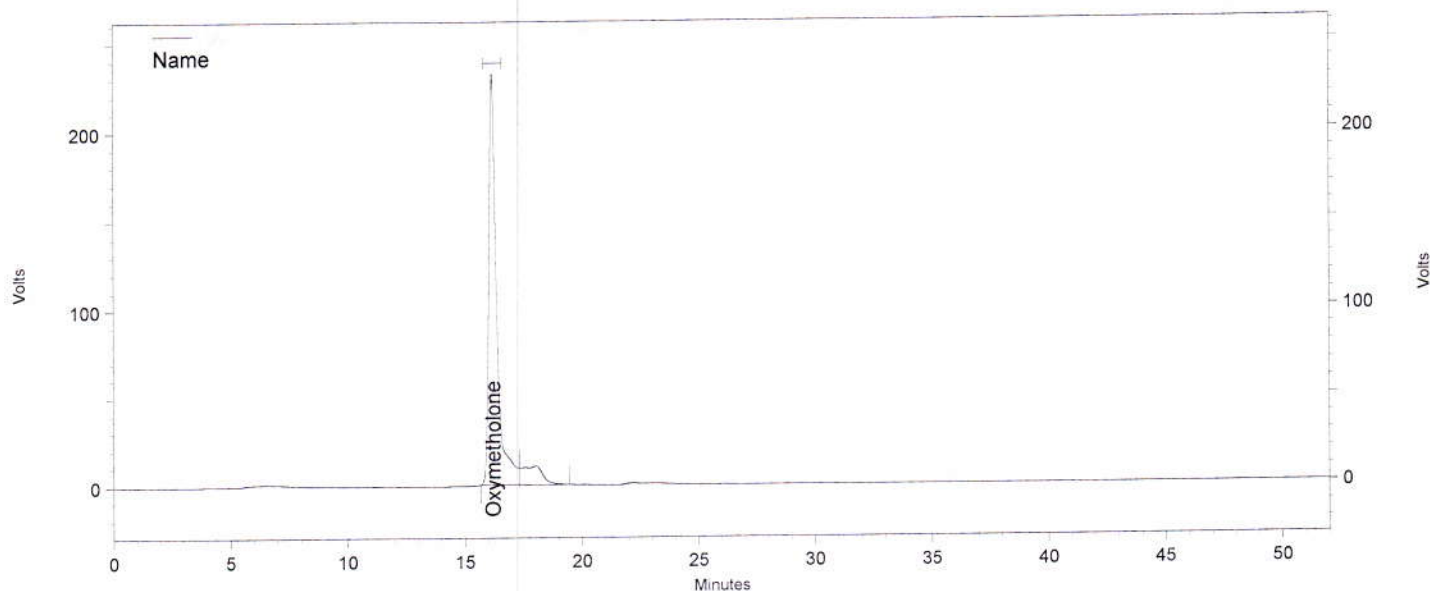

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	16.183	21434648	0.000	89.345
	18.017	2556105	0.000	10.655


Totals		23990753	0.000	100.000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\06.10.2015 02-25-43_C35895 B-Rep1.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
Analysendatum: 06.10.2015 03:20:20
Vial: 4
Probenidentifikation: C34811 B
Operator: admin
Chrom.-Beding. ACN / Wasser, gradient
Säule Nr. 88, 100-5 C18 250 x 4.6 mm
Detektion: UV 240nm

Druckdatum: 06.10.2015 11:46:31
Injektvolumen: 20

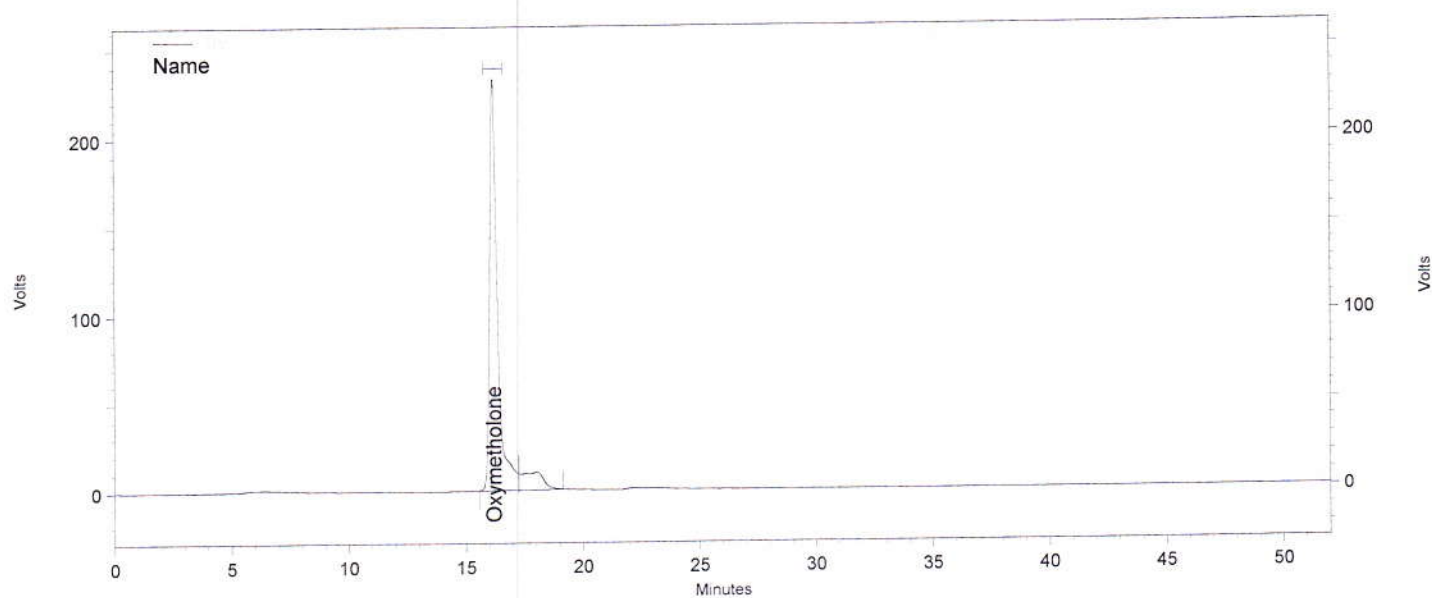

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	16.183	21437295	0.000	89.072
	18.030	2629991	0.000	10.928


Totals		24067286	0.000	100.000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
Steroide\05.10.15\06.10.2015 03-19-07_C35895 B-Rep2.dat
Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met

Report steroid analysis

HPLC: Nr. 6
 Analysendatum: 06.10.2015 04:13:49
 Vial: 4
 Probenidentifikation: C34811 B
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Säule Nr. 88, 100-5 C18 250 x 4.6 mm
 Detektion: UV 240nm

Druckdatum: 06.10.2015 11:46:52
 Injektvolumen: 20

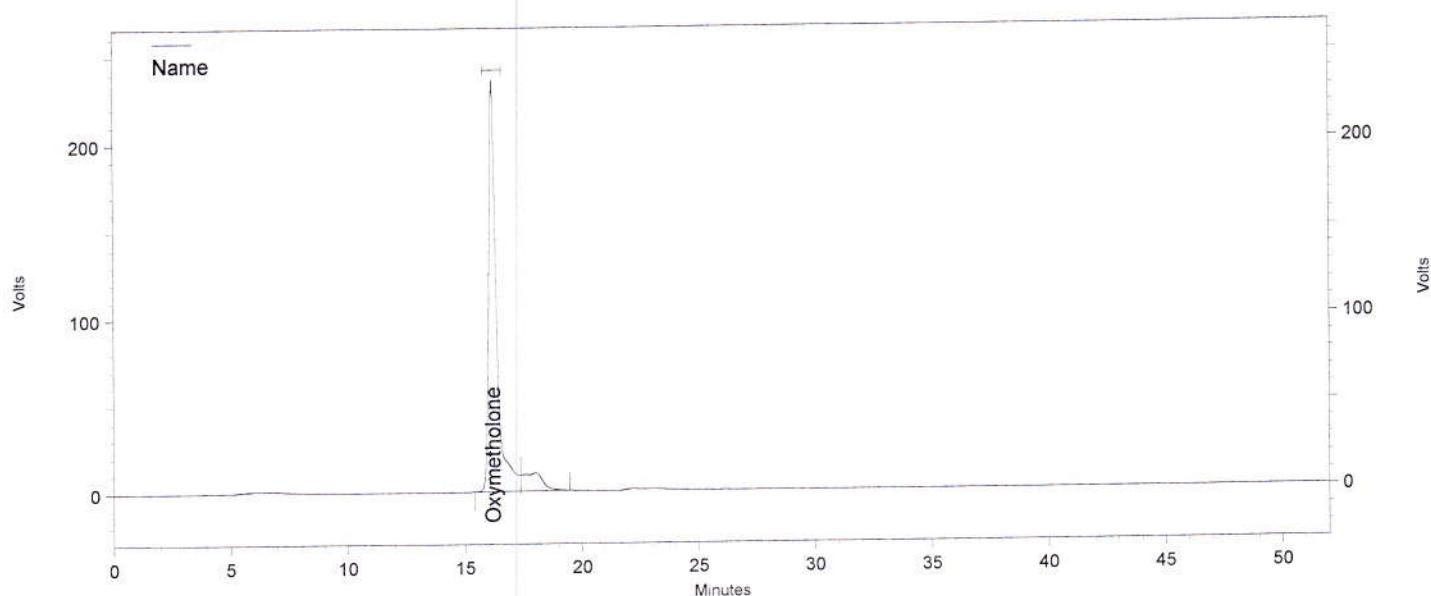

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	16.190	21812112	0.000	89.877
	18.027	2456639	0.000	10.123

Totals		24268751	0.000	100.000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2015\Oktober 2015\Diverse
 Steroide\05.10.15\06.10.2015 04-12-36_C35895 B-Rep3.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_290nm.met